

Reload System of HIFU


# STORMING

## HIFU ANTI AGEING & TIGHTENING



**Non-invasive**  
**No damage on skin surface**  
**Various treatment depth**  
**Minimal risk of adverse events**  
**Natural results- Safe- Effective**  
**No recovery time required**  
**Clinically proven**

**Reduces wrinkles and tightens saggy skin.**  
**Non invasive face lift for cheeks, eyebrows, and eyelids**  
**Redefines jawline and tightens the décolletage.**  
**Natural looking and long-lasting results.**  
**No downtime, safe and effective.**

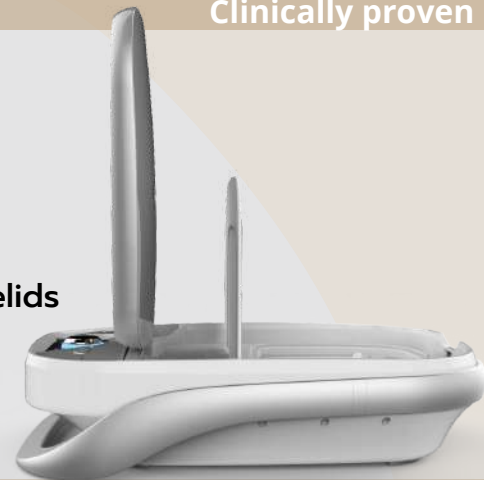


**STORMING High Intensity Focused Ultrasound System**  
**is an innovative non-invasive method for facelifting and skin tightening.**  
**STORMING HIFU technology uses high-intensity focused ultrasound to induce thermal**  
**conduction that shrinks the skin and the subcutaneous layer of fat**  
**metabolism and blood circulation stimulation**  
**new collagen the protein that provides firmness and skin elasticity. The skin is**  
**rejuvenated, looks younger and brighter and the contour is improved**

**Ultrasound can focus energy on target part without harm of surrounding tissue**  
**because it easily passes through the body tissue.**

A safe, effective, non invasive procedure for Face and Body treatments to improve wrinkles and skin laxity of several areas of the face, to reduce the signs of ageing and fat deposits and improve the skin tone in a single session, with potentially long-lasting effect.

**FACE**

**BODY**

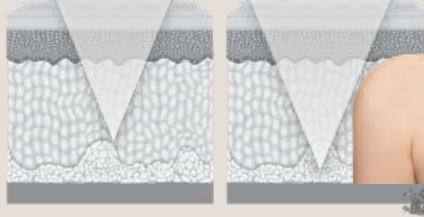

**Upper dermis**

- Forehead
- Periorbital

**Deep dermis**

- Forehead
- Eye area
- Cheek

**SMAS**

- Cheek
- Chin
- Neck

**FAT**

- Abdomen, Plank
- Thigh
- Arm

**FAT**

- Abdomen, Plank
- Thigh
- Arm



**Continuous Working Mode Technology (CWM)**

**Reloads line shots Continuously**

The newest technology in High Intensity Focused Ultrasound. Our 'Storming' HIFU Device can now work continuously. Traditional Cartridges shoot lines one by one, each time starting from the beginning to the end of the cartridge then beginning to end every line shot. The 'Storming' Device keeps shooting continuous shots from left to right and right to left non stop. This makes our device faster with continuous peak power output the whole time non stop.



Portable and

**180° rotary touch panel**

- Clockwise turning up to 180 degrees
- Lightweight and easy to store
- Neat handpiece arrangement
- Display all cartridge information
- Easy to check current shot count



After a STORMING High Intensive Focused Ultrasound treatment

some results will be seen immediately. Most improvement can be seen within 3-4 weeks, with further gains for up to six months.



**Specifications**

Handpiece	1 Handpiece ( Face & Body )
Energy	1.5J MAX ( 1~10 Level )
Spacing	1.0 / 1.5 / 2.0 / 3.0 (mm)
Cartridge Information	4.3Mhz - 3.0mm / 4.5mm (Basic) 4.3Mhz - 1.5mm/ 7mm / 13mm (Optional)
Warranty Shots	Face & Body 10,000 Shots
Voltage / Frequency	100 ~ 240 V / 60 Hz
Dimensions / Weight	410 x 290 x 90 (mm) / 4.2 Kg